



Makers Making Change

A Neil Squire Program



LipSync Wireless

STARTUP GUIDE

Table of Contents

About the LipSync Wireless	3
What Do You Get?	3
Types of LipSync	3
Getting Started	4
Mounting	4
USB Mode and Wireless Mode	5
Mobile Phone and Tablet Compatibility	5
External Power Bank	5
Powering the LipSync	6
Connecting the Device in Wireless Mode	6
Connecting the Device in USB Mode	8
Switching from Wireless mode to USB mode	9
Switching from USB mode to Wireless mode	12
Supporting Applications	13
Using the LipSync	14
Functions	14
Cursor Speed Control	15
Feedback	16
Troubleshooting	16
Joystick Calibration	17

PROGRAM OF:

FUNDING PROVIDED BY:



Foundation



About the LipSync Wireless

The LipSync Wireless is a mouth operated joystick that allows a person who has limited or no control of their hands and arms to control a cursor on a mobile phone, tablet or computer using only head movement. It is an open source assistive technology design by Makers Making Change, a program of the Neil Squire Society that connects volunteer makers with people with disabilities to create more affordable assistive technology.


Open source designs are freely released by their creators for people to make for themselves, modify or incorporate in their own designs for free. Join the Makers Making Change to help people in your community: <http://www.makersmakingchange.com>

Need a build – join us.

Want to make a build – join us.

Want to fund a build for a person with a disability – join us.


What Do You Get?

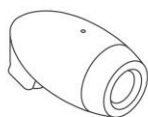


LipSync




Filter x 3

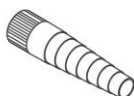

Mouthpiece x 3

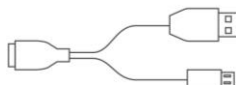

USB adapter

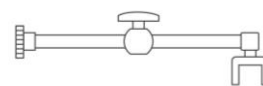

Mounting

Types of LipSync

	Description	Connection Type
LipSync	Works like a mouse	USB
LipSync Wireless	Works like a mouse	USB & wireless
LipSync Macro	Works like a keyboard	USB & wireless
LipSync Gaming	Works like a gaming joystick	USB

1 Getting Started

The LipSync Wireless requires the following items to get started:

- A LipSync Wireless
- Computing device (mobile phone, tablet, or laptop)
- External power bank (optional but recommended)
- Mounting solution for securing the LipSync to a surface

Before plugging the LipSync Wireless, make sure of the following:

1. Computing device is fully charged and/or uses an external power source (e.g. power bank).
2. Computing device and LipSync Wireless are securely mounted to a stable surface.

Mounting

Wheelchair: Every wheelchair has different points available for mounting. Common mounting points available on power wheelchairs are the headrest rail, back of wheelchair, armrest rail, under chair cushion, chair rail, and footplate. Try out different areas and see what works best. If you are not familiar with mounts, consult with a local assistive technology professional.

Desk: Find a clear space on the front edge of your desk or table. Orient the mount bracket toward the user being careful not to move the LipSync. Secure the bracket and adjust the mounting arm so the LipSync is within easy reach. Tighten the knob when the arm is in the desired position.

Checklist for after mounting:

- ☐ Are the clamp and mount secure and do not move?
- ☐ Is there clearance around the wheels, wires, and chair frame with the mount attached?
- ☐ Does the wheelchair fit through doorways?
- ☐ Is the driving view of the user clear with no obstructions?
- ☐ Is the mobile device within viewing angle of the LipSync?
- ☐ Can the user easily reach the LipSync to use it?

USB Mode and Wireless Mode

The LipSync Wireless has two built in modes: USB mouse mode and Wireless mouse mode. By default, the device will be in wireless mode. Users who plan to use the LipSync with multiple electronic devices may benefit from using it in wireless mode whereas users only using it with one device can set the Lipsync to USB mode.

Mobile Phone and Tablet Compatibility

	Operating System	Compatibility
Phone / Tablet	Android ¹	✓
	iOS ²	✓
	Windows	✓
Desktop	Windows	✓
	macOS	✓
	Linux	✓

¹ Only with devices with USB OTG HID (Universal Serial Bus On-The-Go Human Interface Device)

² Only for iOS13+ and iPadOS. For older iOS versions, please use the *LipSync Macro*


Not all mobile phones and tablets support external pointing devices like the LipSync or other mouse replacements. Apple iOS phones and tablets with iOS13 and iPadOS or higher, and only some Android mobile phones and tablets support the LipSync. Check with the manufacturer before purchasing. Look to see if your device supports “USB OTG HID” or “Universal Serial Bus On the Go Human Interface Device” profile.

External Power Bank

The LipSync can be powered using an external power supply like power bank of at least 2700mAH or higher. This prevents the LipSync from drawing battery from your mobile or tablet. The power bank is only supplying power to the LipSync and the mobile or tablet will not be charged.

Powering the LipSync

The order in which you connect the USB cables may affect how the LipSync functions. Plug the cables in the following order:

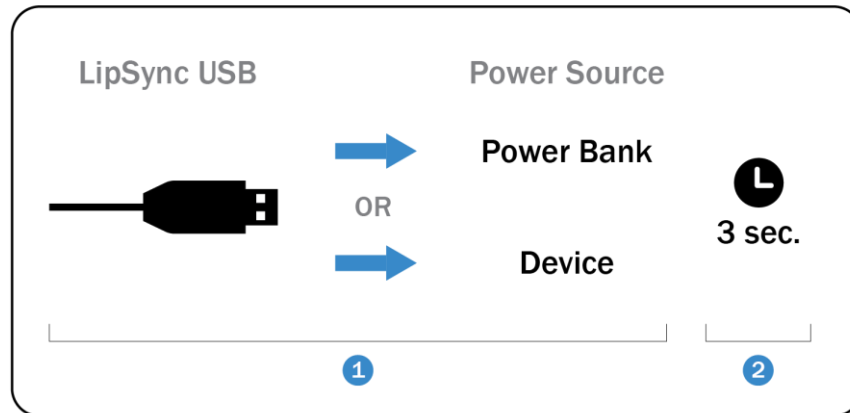

1. Connect the LipSync USB to either a power bank or your device. Connecting to your device for power will draw will cause the device's power to drain quicker than normal.
2. Wait 3 seconds for the LipSync to initialize. Initialization is complete when the LED blinks red and green four times.

Connecting the Device in Wireless Mode

Plugging into a Computer, Laptop or Android & Windows Device

1. Open your device's [Settings](#) > [Bluetooth](#) > Turn it on.
2. Tap on [LipSyncWireless](#) to connect.
3. A [Connected](#) message will show beside the device name when it successfully connects.

Plugging into an iOS Device


1. Open your device's [Settings](#) > [Accessibility](#)




2. Tap on *Touch*





3. Turn on *AssistiveTouch*

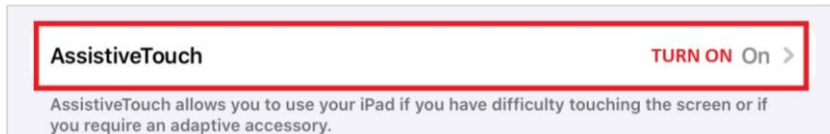

4. Turn off *Always Show Menu*

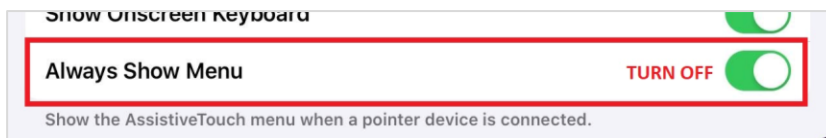

5. Go to *Bluetooth* and turn it on

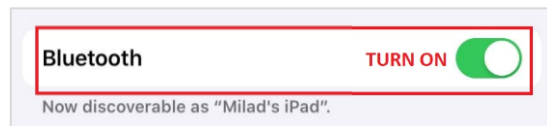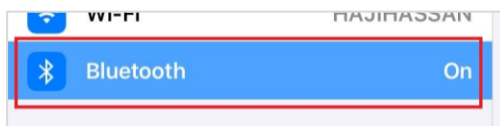


6. Tap on *LipSyncMouse* to connect.



Connecting the Device in USB Mode

Plugging into a Computer, Laptop or Android & Windows Device


The order in which you connect the USB cables may affect how the LipSync functions. Plug the cables in the following order:

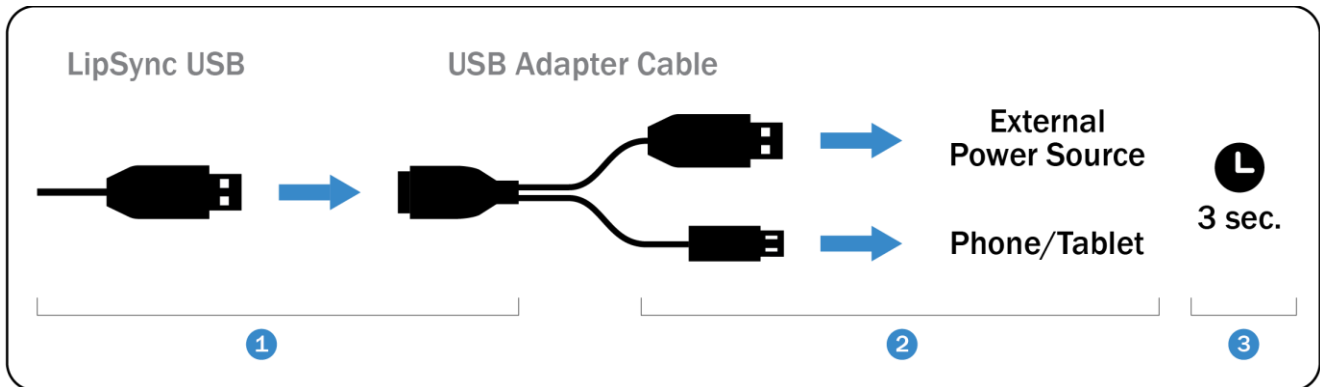

1. Connect the LipSync USB to the USB adapter cable's receptacle.
2. Connect the other end of the adapter into the computing device and external power source.
3. Wait 3 seconds for the LipSync to initialize. Initialization is complete when the LED blinks red and green four times.

Note: If you are using a desktop or laptop the LipSync USB can be directly plugged into the computer USB port.

Plugging into an iOS Device

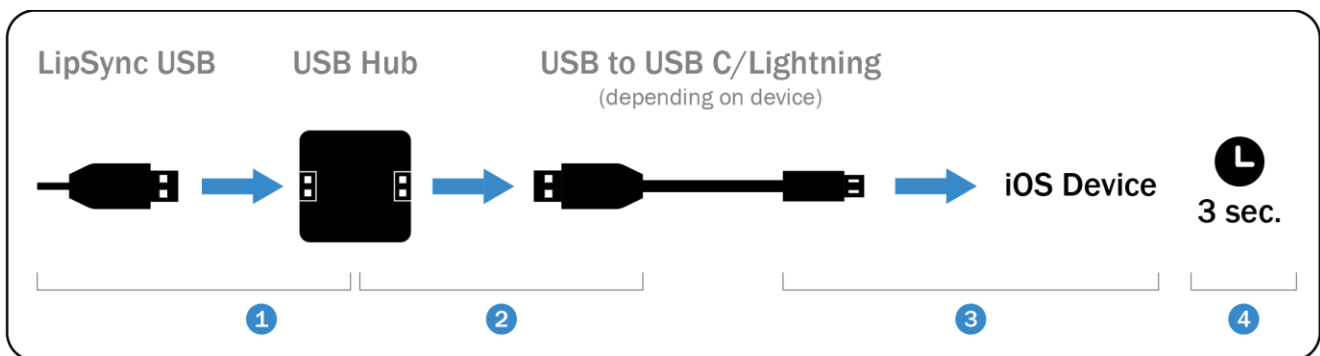

1. Connect the LipSync USB to a USB hub's receptacle.
2. Connect the USB end of a USB to USB C/Lightning cable into the hub's USB receptacle.
3. Connect the other end of the USB to USB C/Lightning cable into the iOS device.
4. Wait 3 seconds for the LipSync to initialize. Initialization is complete when the LED blinks red and green four times.

Switching from Wireless mode to USB mode

1. Remove the threaded adapter and insert.



2. Remove the mouthpiece and filter.




3. Twist off the front LipSync shell.





4. Use your finger and push the internal electronics out from the threaded adapter hole.




5. Gently slide out the internal electronics by wiggling it out.

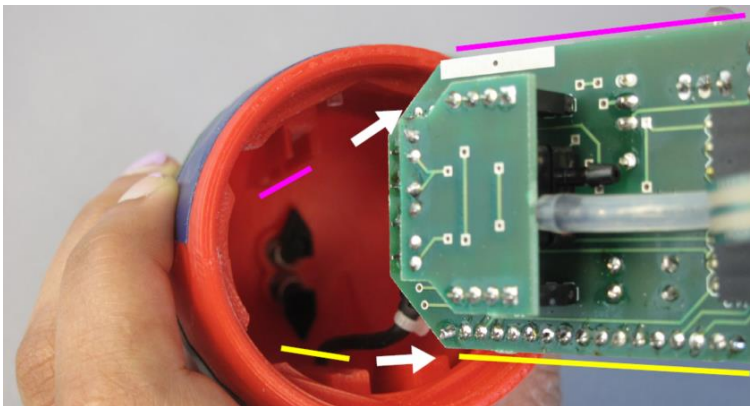

6. Take off the shunt connector and place it sideways.

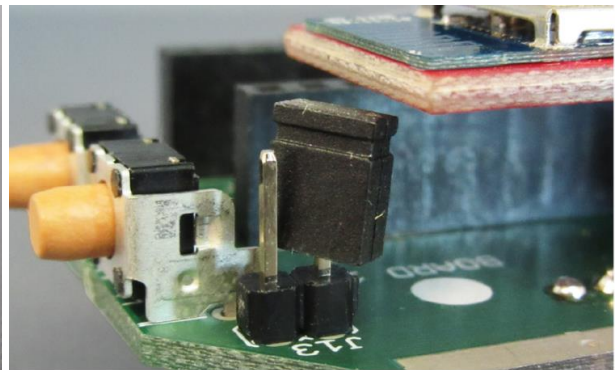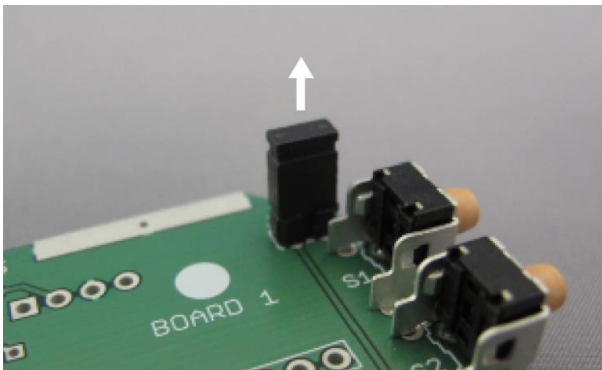

7. Gently insert the internal electronics by wiggling it back into the case. Follow the internal tracks on the inside of the LipSync shell.

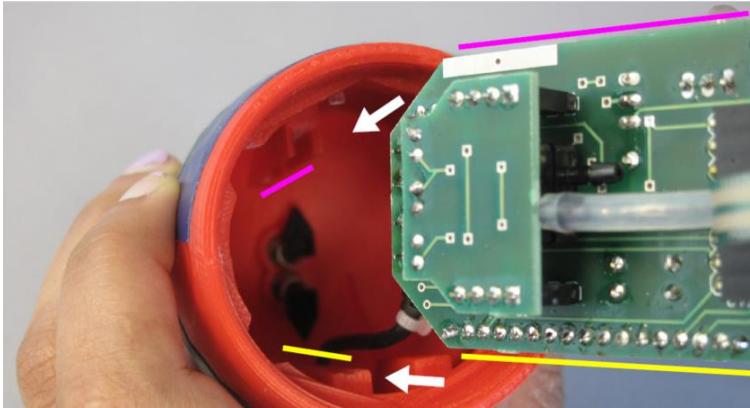

8. Twist on the front LipSync shell



9. Twist on the threaded adapter and insert





10. Twist on the mouthpiece and filter.



Switching from USB mode to Wireless mode


1. Follow steps 1-5 as outlined above.
2. Place the shunt connector onto the two metal leads

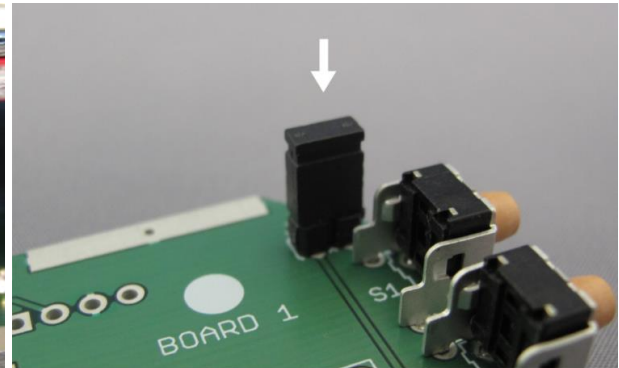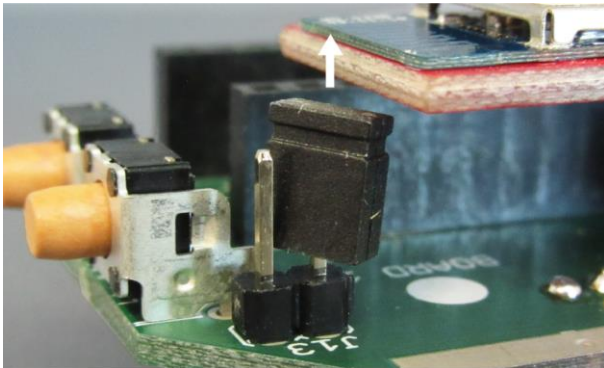

3. Follow Steps 7-10 as outlined above.

Supporting Applications

Some Samsung mobile devices and tablets have an application called **Assistant Menu**, which puts a small menu on the screen that allows quick access to important functions that cannot be accessed without physically touching the device, such as the **Home** button, the **Back** button, **Zoom In and Out**, **Screen Capture**, **Volume**, and **Restart**.

To turn on the Assistant Menu go the **SETTING** application → **MY DEVICE** → **ACCESSIBILITY** → **DEXTERITY AND INTERACTION** → **ASSISTANT MENU**. A little square icon will appear . The icon will always appear on the screen and can be moved by clicking and dragging on it.

If your phone doesn't have **Assistant Menu**, apps like **Assistive Touch For Android**, **Back Button**, **Handy Soft Keys**, **Simple Control** or **Floating Toucher** can be downloaded.

A swiping action is required on some mobile device to unlock it. This can be avoided by setting the phone so that the **SCREEN LOCK** in the Settings Menu is set to **NONE**.

2 Using the LipSync



IMPORTANT: Allow 3 seconds for the LED initialization sequence 🚫 🟢 🚫 🟢 🚫 🟢 to flash **each time** the LipSync is connected to a device before moving the mouthpiece. Moving it before the 3 seconds is complete will affect the usage.

Functions


FOR MOBILE DEVICES AND TABLETS		
Function	Action	Time (seconds)
Tap	Puff	1
Back button ¹	Sip	1
Drag ²	Puff	3
Scrolling ³	Sip	3
Cursor drift reset	Puff	5
Home button ¹	Sip	5

¹ The back/home functions may not work if the device uses on-screen icons. In this case, move the cursor over the icon to tap.

² The red LED will stay on at all times in this mode


³ The green LED will stay on at all times in this mode

Useful function examples:

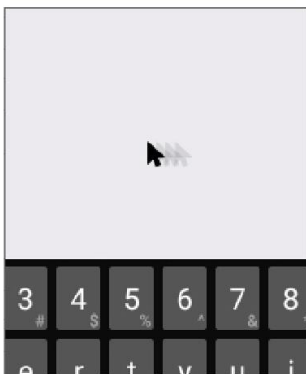


Drift reset / PUFF 5 SEC

- Resolve cursor drift

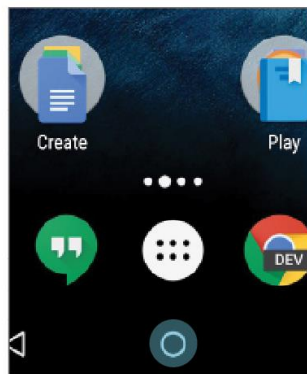


Home button / SIP 5 SEC

- Access home screen

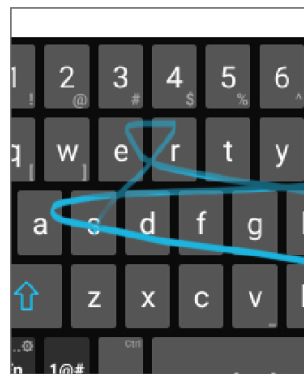

Drag / PUFF 3 SEC

- Unlocking pattern lock
- Keyboard swipe typing

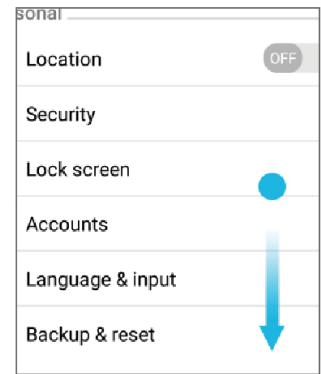


Scrolling / SIP 3 SEC

- Scrollable apps
- Pan through home page

FOR PCs and Laptops		
Function	Action	Time (seconds)
Left click	Puff	1
Right click	Sip	1
Double-click	Puff	1 each
Drag	Puff	3
Scrolling (↑ ↓)	Sip	3
Cursor drift reset	Puff	5
Scrolling (↑ ↓ ← →)	Sip	5

Useful function examples:

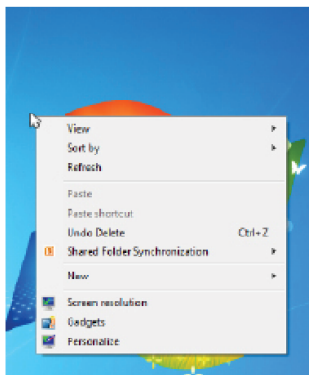


Right click / SIP 1 SEC

- Secondary menu for cut, paste, select all, etc.




Drag / PUFF 3 SEC

- Rearrange icons
- Highlight text

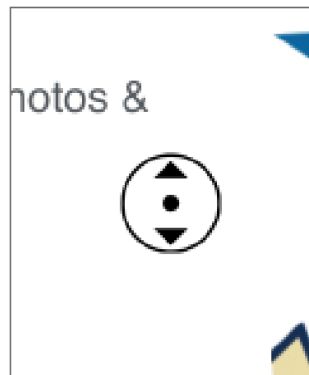

Scroll ↑ ↓ / SIP 3 SEC

- Up and down scroll on a webpage/ document

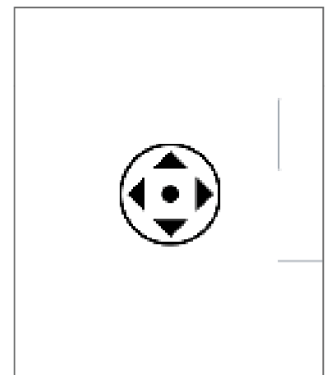

Scroll ↑ ↓ ← → / SIP 5 SEC

- All direction scrolling for large documents

Cursor Speed Control

The LipSync has 11 speed control levels with the default level set to 5 when first initialized.

To adjust the cursor sensitivity and speed, set the computing device's cursor speed to its slowest setting to start, then adjust the LipSync speed control buttons until you get a speed that is comfortable for you. Press the up facing triangular button ▲ on the back of the LipSync to increase the speed and the downward facing button ▼ to decrease. If the speed is still not fast enough, gradually increase the speed on the mobile device or computer.

3 Feedback

Your feedback is very important to us. It helps to support the engineers and designers of the LipSync by learning what is working and what can be improved. It will help inform future development of the LipSync. None of your feedback will be shared publicly without your consent.

Our feedback form can be found at the following link: <https://tinyurl.com/yc57vkpr>

Troubleshooting

Q: Why doesn't the cursor appear on the screen when I plug in the LipSync?

Solution 1: Make a 5 second long puff to trigger the home position reset.

Solution 2: Wait 3 seconds for the LipSync to initialize. Moving the mouthpiece before the 3 seconds is complete will affect the calibration.

Solution 3: If the cursor does not appear after 3 seconds, disconnect the cable from the computing device, reconnect again, and let the LipSync initialize for 3 seconds.

Solution 4: Perform the joystick calibration found on the follow page.

Q: Why does my cursor drift across the screen when I am not moving the LipSync?


Solution 1: Make a 5 second long puff to trigger the home position reset.


Solution 2: Wait 3 seconds for the LipSync to initialize. Moving the mouthpiece before the 3 seconds is complete will affect the calibration.


Solution 3: If the cursor does not appear after 3 seconds, disconnect the cable from the computing device, reconnect again, and let the LipSync initialize for 3 seconds.


Solution 4: Perform the joystick calibration found on the follow page.


Joystick Calibration



1



x 5


x 5


3 sec.
pause

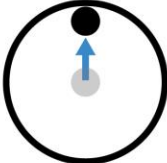
x 6


3


x 1

4

Press both speed adjustment buttons on the back of the LipSync simultaneously for a few seconds

2


LED slowly blinks red & green 5 times


3


When LED blinks green 6 times, move the mouthpiece to the 12 o'clock position (up direction)

4

LED blinks red 1 time when position is measured


x 6

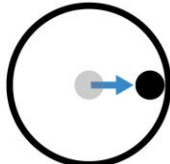
5



x 1


6


When LED blinks green 6 times, move the mouthpiece to the 3 o'clock position (right direction)

6

LED blinks red 1 time when position is measured


x 6

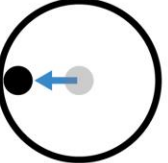
7



x 1


7


When LED blinks green 6 times, move the mouthpiece to the 6 o'clock position (down direction)

8

LED blinks red 1 time when position is measured


x 6

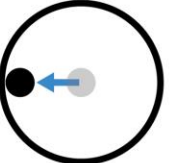
9



x 1


9

When LED blinks green 6 times, move the mouthpiece to the 9 o'clock position (left direction)

10

LED blinks red 1 time when position is measured


x 6


x 6

11

LED will rapidly blink red and green 6 times when calibration is complete